

B

Breast upliftment

V

Vaginal Rejuvenation

P


Pain Management

A

Aesthetics

C

Body Contouring

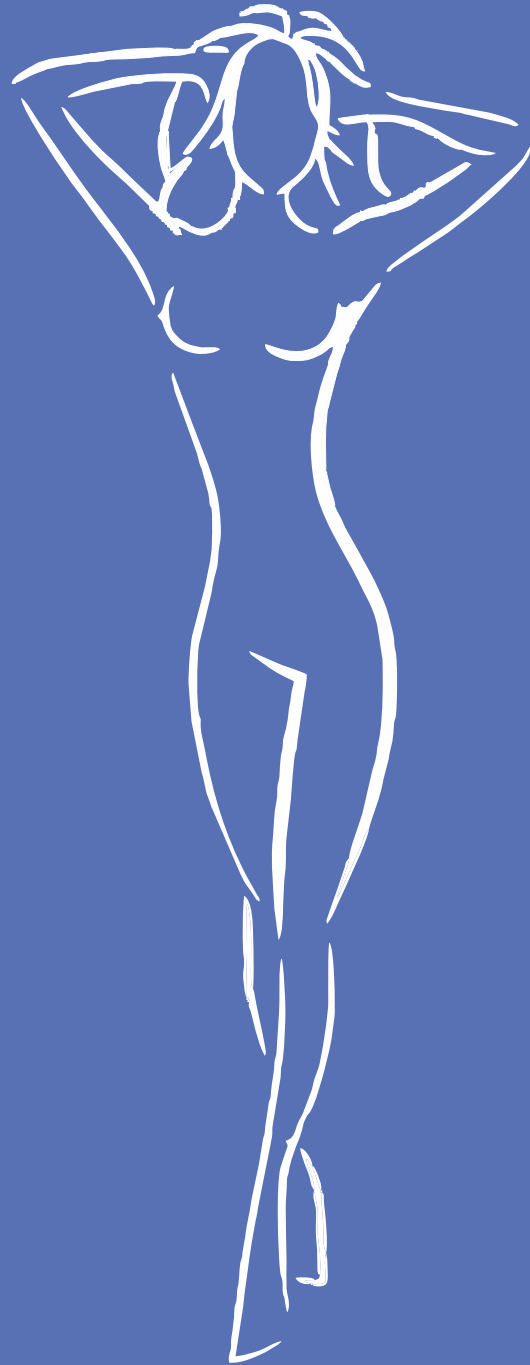

V-juvenation RF device works at a cellular level to stimulate, rebuild and restore collagen within the vaginal tissue - without surgery, without pain and without down time. What makes this treatment unique is that it gently heats the deeper tissue, without causing discomfort to the epidermis.



V - Juvenation
RF - Vaginal
Rejuvenation System

V - Juvenation
for Vaginal Tightening & Urinary Incontinence

3D RF Technology
with Impedance
Feedback and
Temperature Control

The diagram illustrates the V-juvenation RF device, a black handheld wand with a curved tip. Red concentric arcs emanate from the tip, representing the radiofrequency energy being applied. To the left, three circular cross-sections of the vaginal canal are shown. The top one is labeled 'BEFORE' and shows a wider, less defined canal. The middle one is labeled 'DURING' and shows the device's tip inserted into the canal. The bottom one is labeled 'AFTER' and shows a narrower, more defined canal. To the right, a list of benefits is provided.


Benefits

- Vaginal (Internal) Tightening
- Vulvar/Labial (External) Tightening
- Decreased Vaginal Dryness
- Reduced Urinary Incontinence
- Increased Sexual Satisfaction

Aesthetics: RF energy is an excellent remedy against wrinkles and skin dehydration. By stimulating cellular capabilities, it reactivates mature skin,, restoring elasticity and softness




Pain Management Divlabs PAIN FREE Radiofrequency Therapy distributes high frequency electromagnetic energy into the desired tissue which leads to selective tissue hyperthermia. Thanks to this mechanism patients experience instant pain relief, muscle relaxation, oedema reduction, tissue regeneration support and healing.

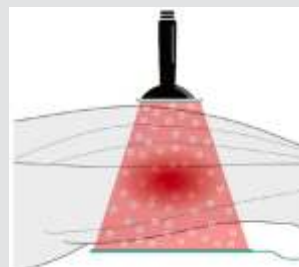

Body Contouring & Sculpting



3D RF increases the blood supply from sub dermal to dermal region (muscles to fat)

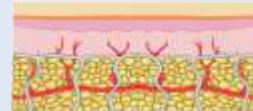

Enhanced blood circulation burns the target fat equivalent to an extensive exercise



3D RF HEXA
(For Abdomen, Hips, Buttocks and Thighs)



3D RF QUAD
(For Calves, Arms, Neck and Shoulder)



3D RF TRIAD
(For Aesthetics Face, Chin and Localised Fat)

- Vivid Shape uses a combination of radio frequency energy (RF) infrared (IR) technology and mechanical tissue manipulation for effective fat burning instantaneously.
- The hot spot created by the 3D RF deep heats the fat cells, connective tissue & dermal collagen fibers resulting in an increased rate of cellular metabolism and facilitating lymphatic drainage.
- Neo Collagenesis - The combination of 3D-RF and EMC increases the collagen production resulting in smoothing of skin texture and reduction in skin laxity and wrinkles.

VIVID SHAPE

General

Detection Control System
Auto Temperature Control
Display Touch Screen 10,4"

Power source	110 - 240 VAC 50-60 Hz
Device type	BF
Insulation class	I
Medical Device	Class IIb - CE
Casing	Fire proof ABS
Fuses	2 x 250 V T4A
RF output power	50 W (aesthetics)
RF output power	300 W (medical)
Output frequency	1.2 MHz.
Connectivity	Wifi, LAN, USB
Size l x h x d	400 x 310 x 230 mm
Accessories	Cart, trolley
Weight	10 kg
Accuracy of displayed values	± 20%

DIVLABS HEALTHCARE PVT. LTD.

Regd. Off.: Plot No. 17, 1st Floor, Kh. No. 44, Toganpura, Sector-9, Dwarka, New Delhi-110077, INDIA
Office: 320, City Centre Mall, Sector-12 Metro Station, Dwarka, New Delhi-110078, INDIA

Tel. : 91-11-28081996 | Mobile : **+91 9811259555, 9811222377**
E-id: ceo.divlabs@gmail.com | Web: www.divlabsindia.com